


Vented I.V. Giving Set for Paclitaxel

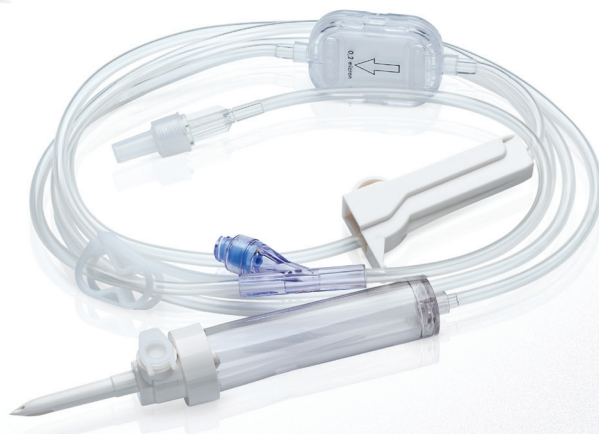

Paclitaxel Family of Sets

The Paclitaxel sets are infusion sets with a 0.2 micron filter for the administration of cytotoxic drugs requiring an in-line filter.

Vented I.V. Giving Sets. On the right, with Luer Activated Injection Site.

TECHNICAL DATASHEET

Vented I.V. Administration (Giving) Set for Paclitaxel

The vented IV administration sets offer customized sets for administration of Paclitaxel or other viscous drugs or drugs that require in-line filtering. The vented IV administration set for Paclitaxel with 0.2 micron filter and vented spike allows comfortable and simple Paclitaxel administration via a luer lock connection. This set replaces the need for an external filter and a separate IV administration set. The Paclitaxel general infusion sets provide luer lock access points to the patient's IV port. The Paclitaxel sets have an option for an additional needleless injection site.

For a closed drug administration, connect a Tevadaptor® Syringe Adaptor/ Syringe Adaptor Lock to the end of a vented IV administration set and a Tevadaptor® Luer Lock Adaptor to the patient's port.

TECHNICAL DATASHEET

VENTED I.V. GIVING SET FOR PACLITAXEL

VENTED I.V. GIVING SET FOR PACLITAXEL WITH LUER ACTIVATED INJECTION SITE

Order Information

Article number / Reference	MG245537	MG245538
User box quantity	75 pcs	75 pcs
Transport box quantity	75 pcs	75 pcs

Product Characteristics

Product Technical data		
Tubing length	180 cm	180 cm
Volume of the tubing	19 ml (approx.)	19 ml (approx.)
Spike type	Vented 2-way spike	Vented 2-way spike
External filter size	0.2 micron	0.2 micron
Proximal end of the tubing	Male luer lock	Male luer lock
Clamp type	Pinch clamp and roller clamp	Pinch clamp and roller clamp
Injection site sealing	Regular injection site	Luer activated injection site (Needle-less)
Shelf life	5 years	5 years
Sterilization	Ethylene Oxide	Ethylene Oxide
Drip chamber and drip tube	1ml=20 drops of distilled water	1ml=20 drops of distilled water

Material Specifications

Component	Material	
Protective cap on spike	Polyethylene	Polyethylene
2-way Spike	ABS	ABS
Injection site sealing	Polyisoprene and ABS	Polycarbonate and silicone
Soft tubing	DEHP-free PVC	DEHP-free PVC
Pinch clamp on tubing	Acetal	Acetal
Male luer lock body	ABS	ABS
Drip chamber	DEHP-free PVC and ABS body Polyamide 6.6 filter	DEHP-free PVC and ABS body Polyamide 6.6 filter
Roller clamp	ABS	ABS
0.2 Micron filter	MBS-Cyrolite, PES- Polyethersulfone and PTFE	MBS-Cyrolite, PES- Polyethersulfone and PTFE

DEHP-free PVC. Non-pyrogenic. Free from latex, lead, cadmium, arsenic and chrome. Meets standards and regulations with regards to the environment.

Regulation

CE 0483
Disposal of components should be according to local regulations and procedures for handling of hazardous drugs

Symbols taken from ISO 15223-1:2016



Sterilized using ethylene oxide

Non-pyrogenic

Do not re-sterilize

Do not re-use

Do not use if package is damaged



Keep dry



Keep away from sunlight



Donawa Lifescience Consulting Srl
Piazza Albania 10,
00153 Rome, Italy

CE 0483



Simplivia Healthcare Ltd.
North Industrial Zone
P.O. Box 888
Kiryat Shmona, Israel 1101801



www.simplivia.com



info@simplivia.com



Follow us

