



# Syringe Adaptor Sterile Replacement Cap

## TECHNICAL DATASHEET

The Chemfort™ Syringe Adaptor Sterile cap is a PVC-free component and fits Chemfort™ Syringe Adaptor and Chemfort™ Syringe Adaptor Lock Elastomer side.


The sterile protective cap included in the package, can be reused for capping after activation. Note that the cap must be wiped with 70% IPA or alcohol before use.

**Notice** the cap is not sterilized after opening.


The recapping action allows safe delivery of a drug-filled syringe to the ward in case of bolus injections, or while administering through IV set or Chemfort™ Closed Administration Set in the tail end.

### Indications for Use:

Chemfort™ is a Closed System Transfer Device (CSTD) that mechanically prohibits the release of drugs, including antineoplastic and hazardous drugs, in vapor, aerosol, or liquid form during preparation, reconstitution, compounding and administration, minimizing exposure of individuals, healthcare personnel and the environment to hazardous drugs. Chemfort™ prevents the introduction of microbial and airborne contaminants into the drug or fluid path for up to 7 days.

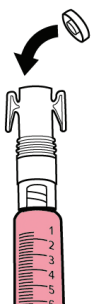

Chemfort™ Syringe Adaptor for bolus injection

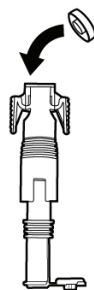

On Chemfort™ Closed Adaptor Spike Port (connected to the IV set spike)



Chemfort™ Closed administration Secondary set



Each cap is packed in a separate blister packaging. Each selling box contains 100 units connected in strips of 5 sterile self-packed caps.

## Syringe Adaptor Sterile Replacement Cap

### TECHNICAL DATASHEET

#### Order information

Article number / Reference	MG245268
User box quantity	100 pcs
Transport box quantity	400 pcs

#### Product Characteristics

Tolerates disinfection with alcohol 70% and chlorhexidine in alcohol (5 mg/70% alcohol)

#### Product Technical data

External fitting	For re-capping Chemfort™ Syring Adaptor and Syringe Adaptor Lock
Cap weight	0.412 gr
Shelf life	3 years
Sterilization	Ethylene Oxide

#### Material Specifications

DEHP-free PVC. Non-pyrogenic. Free from latex, lead, cadmium, arsenic and chrome. Meets standards and regulations with regards to the environment

Component	Material
Syringe Adaptor Replaceable cap	Polyethylene

#### Regulation

CE 0483

FDA 510(k) Cleared

Disposal of components should be according to local regulations and procedures for handling of hazardous drugs

Symbols taken from ISO 15223-1:2021



Consult instructions for use



Sterilized using ethylene oxide



Non-pyrogenic



Do not re-sterilize



Do not re-use



Do not use if package is damaged



Keep dry



Keep away from sunlight



Donawa Lifescience Consulting Srl  
Piazza Albania 10,  
00153 Rome, Italy

CE 0483

SIMPLIVIA™

Simplivia Healthcare Ltd.  
North Industrial Zone  
P.O. Box 888  
Kiryat Shmona, Israel 1101801



[www.simplivia.com](http://www.simplivia.com)



[info@simplivia.com](mailto:info@simplivia.com)



Follow us