


TevaTree Inline Set

4 Female Luer Connectors

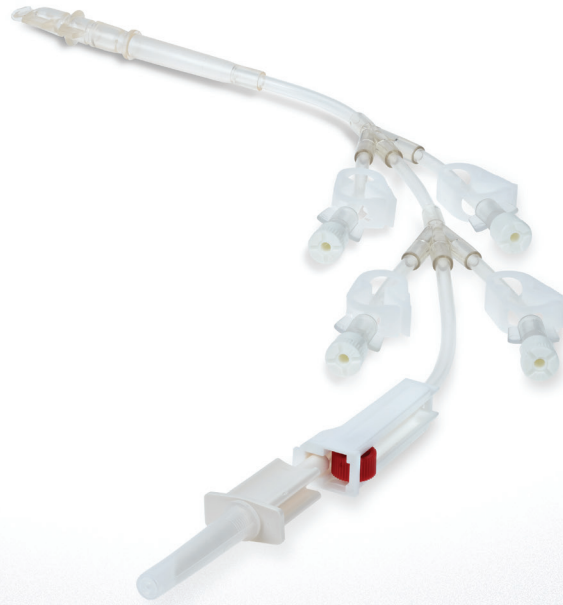

TevaTree
Administration
Family of Sets

TECHNICAL DATASHEET

TevaTree Family

Offers sets with multiple arms for connection to multiple cytotoxic drugs (via a Tevadaptor® Connecting Set) simultaneously with a saline IV bag for either full or in-line administration.

TevaTree Inline Set, 4 Female Luer Connectors:

The TevaTree in-line general set offers an in-line set with 4 arms to connect multiple cytotoxic drug IV bags and a saline IV bag. The TevaTree arms are equipped with female luer lock connectors to safely connect to Tevadaptor® Connecting Sets, minimizing the risk of leakage or disconnection.

The TevaTree inline set can be connected to a general IV set.

When used with Tevadaptor® components, the TevaTree general set is part of a full solution for your administration needs.

TECHNICAL DATASHEET TEVATREE INLINE SET, 4 FEMALE LUER CONNECTORS:

Order Information

Article number / Reference	MG245540
User box quantity	50 pcs
Transport box quantity	50 pcs

Product Characteristics

Product Technical data	
Tubing length	40 cm
Volume of the tubing	4 ml
Spike type	Non-vented spike
Clamp type	Pinch clamp and roller clamp
Shelf life	5 years
Sterilization	Ethylene Oxide
4 female luer connectors	For secondary sets connection
Proximal end of tubing	Twist off connector and internal protective membrane

Material Specifications

DEHP-free PVC. Non-pyrogenic. Free from latex, lead, cadmium, arsenic and chrome. Meets standards and regulations with regards to the environment.

Component	Material
Protective cap on spike	Polyethylene
1-way Spike	ABS
Soft tubing and connectors	DEHP-free PVC
Pinch clamp on tubing	Acetal
Roller clamp	Polystyrol and Polystyrene
Spike port adaptor (twist off connector)	DEHP-free PVC

Regulation

CE 0483
Disposal of components should be according to local regulations and procedures for handling of hazardous drugs

Symbols taken from ISO 15223-1:2016



Sterilized using ethylene oxide



Non-pyrogenic



Do not re-sterilize



Do not re-use



Do not use if package is damaged



Keep dry



Keep away from sunlight



Donawa Lifescience Consulting Srl
Piazza Albania 10,
00153 Rome, Italy

CE 0483

SIMPLIVIA™

Simplivia Healthcare Ltd.
North Industrial Zone
P.O. Box 888
Kiryat Shmona, Israel 1101801



www.simplivia.com



info@simplivia.com



Follow us

