


TevaTree Inline Set

4 Female Luer Connectors

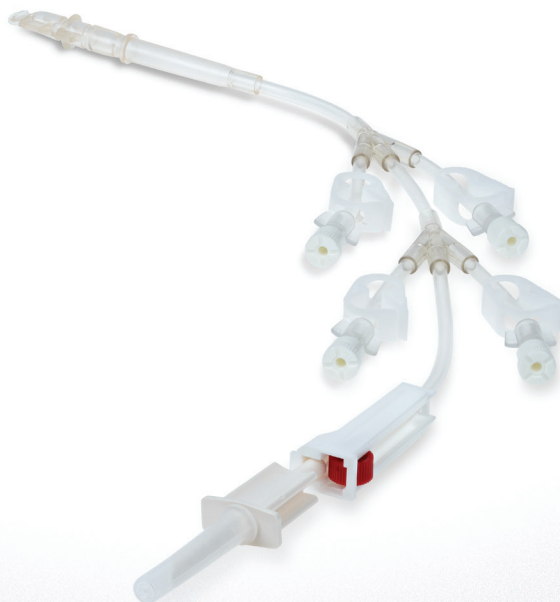

TevaTree
Administration
family of sets

TECHNICAL DATASHEET

TevaTree Family

Offers sets with multiple arms for connection to multiple cytotoxic drugs (via a Tevadaptor® Connecting Set) simultaneously with a saline IV bag for either full or in-line administration.

TevaTree Inline Set, 4 Female Luer Connectors:

The TevaTree in-line general set offers an in-line set with 4 arms to connect multiple cytotoxic drug IV bags and a saline IV bag. The TevaTree arms are equipped with female luer lock connectors to safely connect to Tevadaptor® Connecting Sets, minimizing the risk of leakage or disconnection.

The TevaTree inline set can be connected to a general IV set.

When used with Tevadaptor® components, the TevaTree general set is part of a full solution for your administration needs.

TECHNICAL DATASHEET

TEVATREE INLINE SET, 4 FEMALE LUER CONNECTORS:

Order Information

Article number / Reference	1MG245540
User box quantity	50 pcs
Transport box quantity	50 pcs

Product Characteristics

Product Technical data	
Tubing length	40 cm
Volume of the tubing	4 ml
Spike type	Non-vented spike
Clamp type	Pinch clamp and roller clamp
Shelf life	5 years
Sterilization	Ethylene Oxide
4 female luer connectors	For secondary sets connection
Proximal end of tubing	Standard spike port with twist off connector and internal protective membrane

Material Specifications

Component	Material
Protective cap on spike	Polyethylene
1-way Spike	ABS
Soft tubing and connectors	DEHP-free PVC
Pinch clamp on tubing	Acetal
Roller clamp	Polystyrol and Polystyrene
Spike port adaptor (twist off connector)	DEHP-free PVC

DEHP-free PVC. Non-pyrogenic. Free from latex, lead, cadmium, arsenic and chrome. Meets standards and regulations with regards to the environment.

Regulation

CE 0483



Consult instructions for use

STERILE EO



Non-pyrogenic



Do not re-sterilize



Single use only, do not re-use



Not made with DEHP



Not made with natural rubber latex



Do not use if package is damaged



Keep dry



Keep away from sunlight

EC REP

CE 0483



TEVA Medical Ltd.
North Industrial Zone
P.O. Box 888 Kiryat Shmona,
Israel 1101801

Disposal of components should be according to local regulations and procedures for handling of hazardous drugs.